

AH Plus® Bioceramic Sealer

Step-by-step Guide

- 1 Disinfect and dry the root canal after preparation is completed
- 2 Choose the master cone that fits snugly the apical preparation size
- 3 Remove the grey syringe cap and extrude a small amount of sealer and discard
- 4 Place a new 24-gauge tip on the syringe
- 5 Check if the sealer can be extruded smoothly



Properties, Warnings and Precautions:

- Setting Time: 2-4 hours
- Do not use in patients allergic to any of the components of AH Plus® Bioceramic Sealer.
- Avoid overfilling
- Sealer should be used immediately after placing the 24-gauge tip because it can harden when in contact with air
- Storage: After use clean the outside of the syringe and remove any excess sealer. Twist the storage cap provided with the product tightly onto the syringe head. Place the syringe into the foil pouch and ensure to seal the pouch


Distributed By:

Maillefer Instruments Holding Sàrl, Chemin du Verger 3,
CH-1338 Ballaigues, Switzerland






- 6 Insert the 24-gauge tip no further than the middle third of the root canal
- 7 Inject the sealer until it is visible at the root canal orifice. Maintain the tip immersed in the sealer during injection to minimize the inclusion of voids

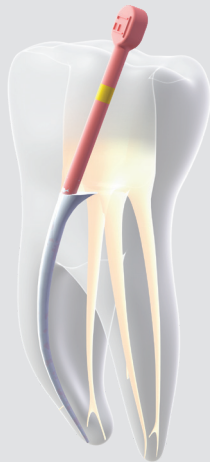

- 8 Insert the master cone into the root canal and push it to the apical stop
- 9 Avoid excessive pressure to minimize or prevent extrusion beyond the apical foramen

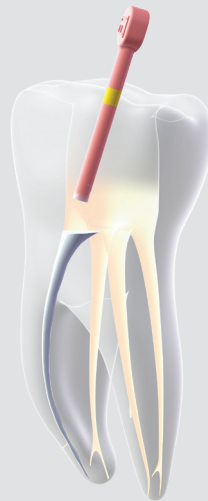

- 10 Cut and remove the coronal portion of the master cone at the root canal orifice
- 11 Compact the coronal portion of the cone with an appropriately sized and fitted plugger
- 12 Verify with radiographs

dentsplysirona.com/ifu

INL / B EN AH00 TIP 000 / 00 / 2021 - created 06 / 2021

EC REP

Biogeri Am Pfeilshof 12, DE-22303 Hamburg, Germany



Manufactured By:

Maruchi, 2-208, Medical Industry Complex Bldg.,
42-10, Taejanggongdangil, Wonju-si, Gangwon-do, KR-26311

